

PHILIPS

DreamStation 2


CPAP Advanced

Auto CPAP Advanced



© 2021 Koninklijke Philips N.V. All rights reserved. Specifications are subject to change without notice. Trademarks are the property of Koninklijke Philips N.V. or their respective owners.

www.philips.com/respironics

 Respironics Inc.
1001 Murry Ridge Lane
Murrysville, PA 15668 USA

REF 1149265
1149265 R02
ZL 02/15/2021

Patient

quick start guide



Connecting to DreamMapper

Sync your device to your tablet or mobile phone using the DreamMapper app, available in the App Store and Google Play.




1

With the device and mobile phone powered up and in close proximity, initiate **Bluetooth** setup and pairing in **DreamMapper**. Follow the pairing instructions in **DreamMapper**.




2

A **6-digit PIN** will appear on the DreamStation 2 screen. Depending on your phone, you will need to either enter the **PIN** into **DreamMapper** or tap **"Pair"** on the device screen within 30 seconds to complete the pairing.

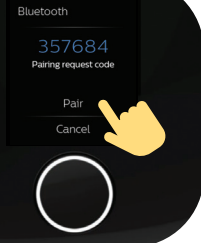


3

The screen displays **"Pairing to Device"** and then **"Success: Device is now paired."**

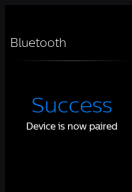

Get started with therapy

To begin therapy, press the **Therapy button**.
A **blue ring** indicates therapy is on.

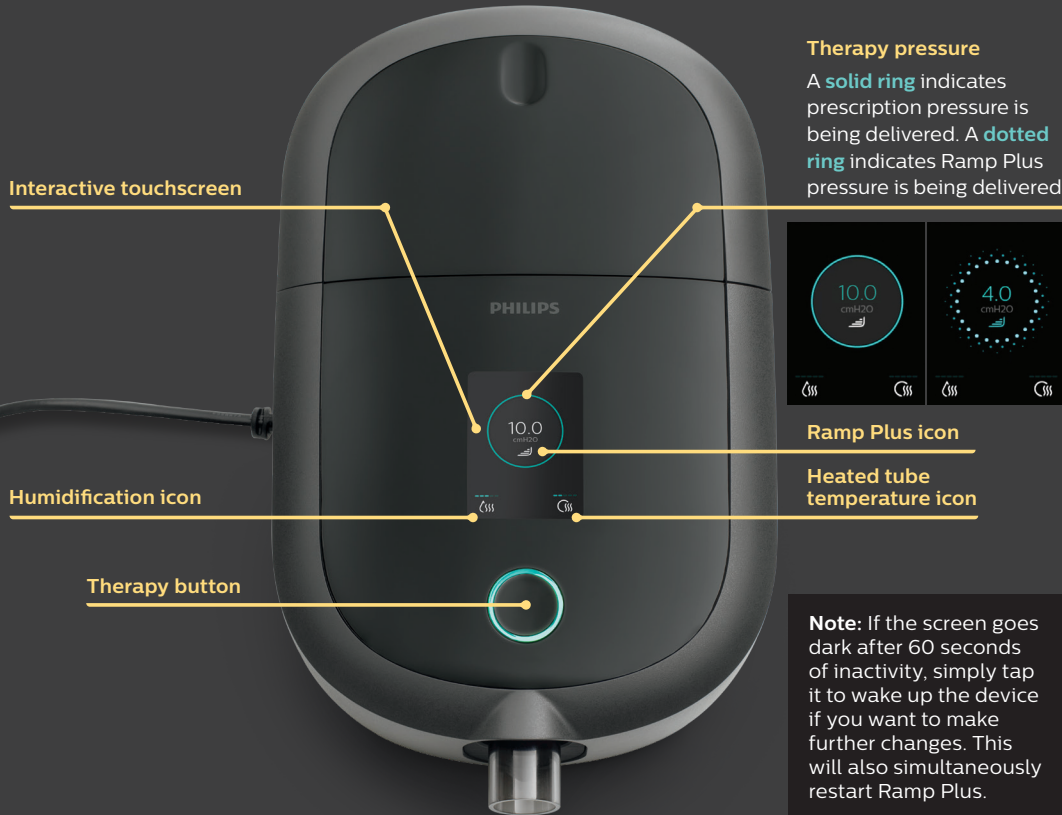

Customize your settings

Ramp Plus

The **Ramp Plus** feature allows you to set a comfortable starting pressure while you fall asleep.

1. Tap the **Ramp Plus** icon to view and adjust your **Ramp Plus** pressure.
2. Tap the **+ or - symbols** until you reach a comfortable setting.
Ramp Plus settings can be turned off or set from 4 to 10 cm H₂O.
3. Tap “**Confirm**” to save the setting.


Once **Ramp Plus** is set, it will be saved and automatically activated for future therapy sessions without tapping the **Ramp Plus** icon again.



Humidification (if available)

1. Tap the **humidification** icon to view your current setting.
2. Tap again until you reach a preferred moisture level.
Humidification can be turned off or set from 1 (minimum) to 5 (maximum).


Once **humidification** is set, it will be saved and automatically activated for future therapy sessions.

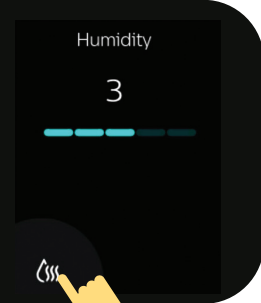

Heated tube temperature (if available)

1. Tap the **heated tube** icon to view your current temperature setting.
2. Tap again until you reach your preferred temperature. Heat can be turned off or set from 1 (minimum) to 5 (maximum).

Once the **temperature** is set, it will be saved and automatically activated for future therapy sessions.

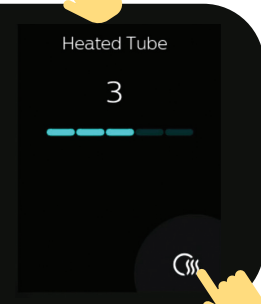


Therapy data summary

To stop airflow, press the **Therapy button**. A summary of your therapy data will be displayed. Access to certain information will depend on what your provider has set.


Last Session screen

Information in **blue** indicates that you can tap to learn more or take action.

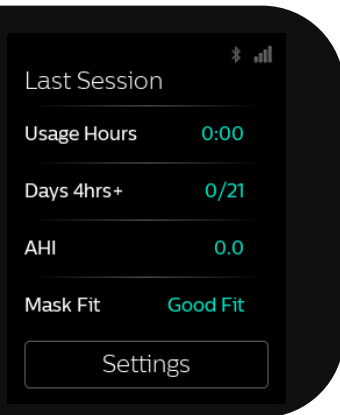


1

The **Usage Hours screen** displays the amount of time you received therapy in hours and minutes during your last 3 sessions.

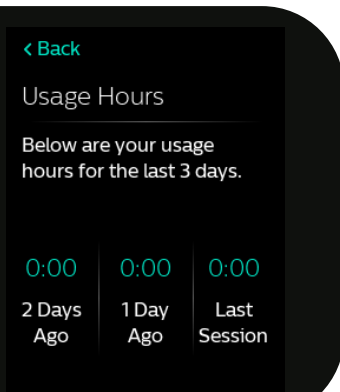


2

The **Days 4hrs+ screen** displays the total number of sessions lasting 4+ hours that you've had in the past 30 days. Once you have reached your goal, "Complete" will be displayed.

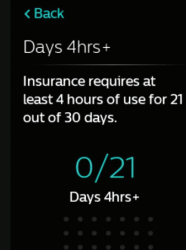

3

The **AHI screen** displays the Apnea-Hypopnea Index (AHI) value for the most recent one-day time frame (if enabled by your provider).

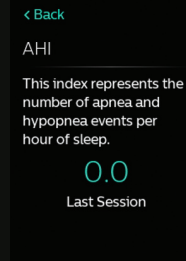

4

The **Mask Fit screen** displays how your mask fit during your last session. You can also run a Mask Fit check. Put on your mask and tap "Start" to initiate the Mask Fit check.




Settings

Comfort, device, and data


Access to certain information will depend on what your provider has set.

- 1** To make adjustments to settings and to access additional information, tap “Settings.”

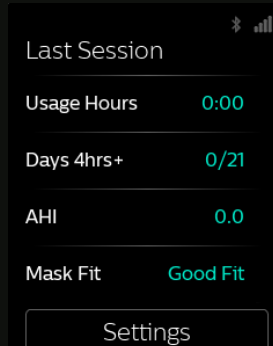

- 2** Adjust available **Comfort settings**. Information in **blue** indicates that you can tap to make changes.


- a. Ramp Plus time and pressure
- b. Humidity
- c. Heated tube (if available)
- d. Tube size
- e. FLEX

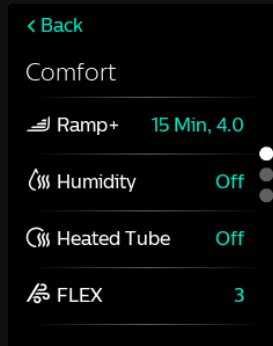

3

Swipe up to view additional device information like adjusting airplane mode, clearing Bluetooth pairings, and resetting time or language.


Tap on the information in **blue** to make changes. You can also view the device's serial number, software version, and modem status.

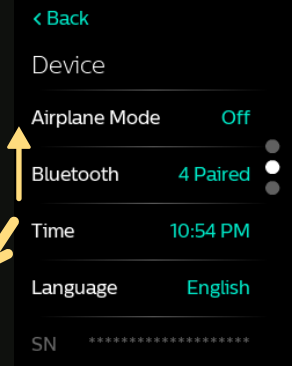

4

Tap “View Reports” to select a date range and view additional therapy data summaries.

To exit out of **Settings** at any time, tap “Back” to return to the **Last Session** screen.

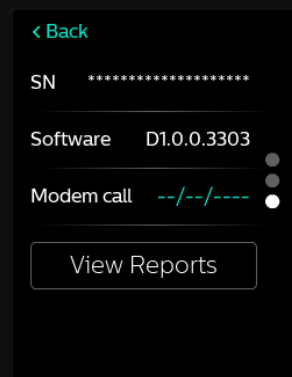

Consult your User Manual for detailed troubleshooting information.

Cleaning and maintenance

To keep your device and components in good working condition, you must clean them regularly. It's also important for maintaining your good health.


Humidifier water tank, seal, and lid

Clean daily with water and mild dishwashing detergent.



Tubing and adaptors

Clean weekly with water and mild dishwashing detergent.

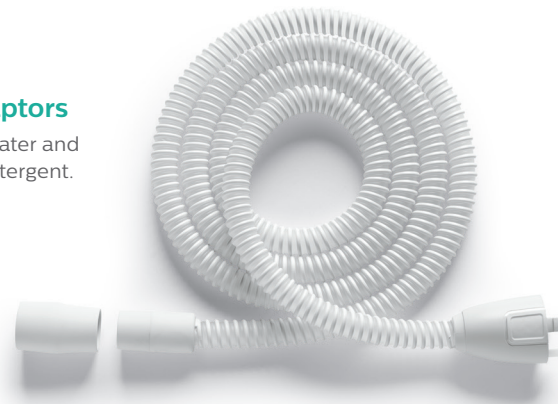

Filters

Clean gray, reusable filter every two weeks with water. Replace every 6 months.

Do not clean blue, disposable filter. Replace every 30 nights.



Consult your User Manual for detailed care and cleaning information.